




AP-Prime II

INSTRUCTION MANUAL

MASSAGE CHAIR

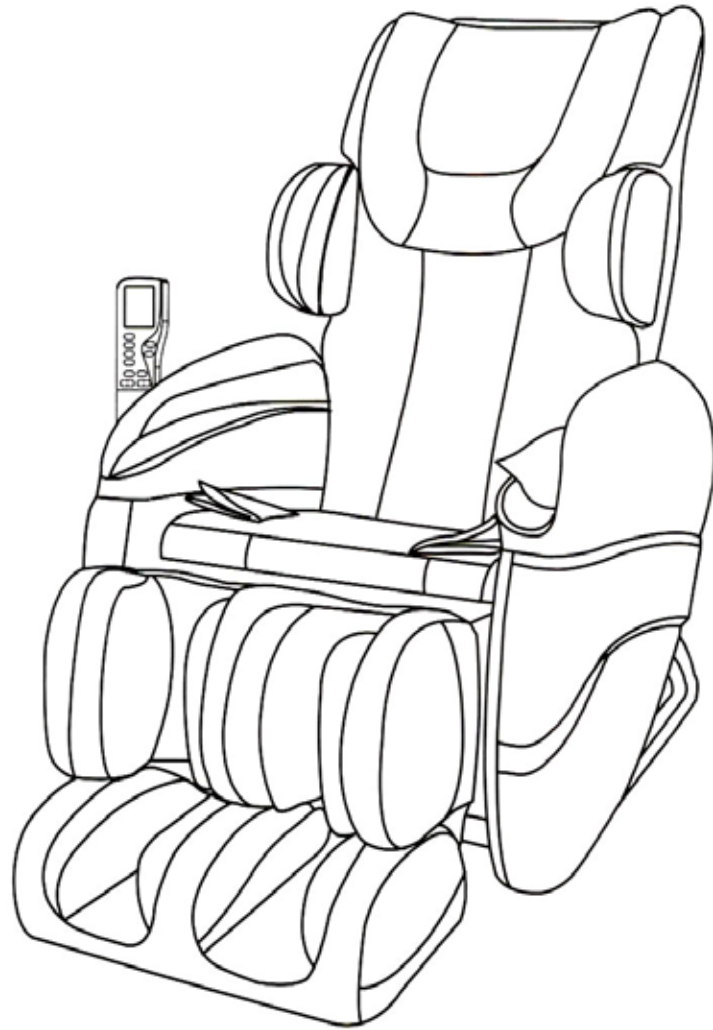

BE SURE TO READ THIS INSTRUCTION MANUAL CAREFULLY BEFORE OPERATING THE CHAIR.
PLEASE KEEP THIS MANUAL IN SAFE KEEPING FOR FUTURE REFERENCE

CONTENTS

Thank you for purchasing our massage chair, hopefully the product will provide will meet your massage needs. For service issues, please contact your local seller or manufacturer, never dismantle or repair this product by an unauthorized technician to avoid any injury or product damage.

I	IMPORTANT SAFETY INSTRUCTIONS	1-3
	PARTS NAME.....	4
I	INSTRUCTIONS ON FUNCTIONS OF THE CONTROL PANEL	5-7
I	HOW TO OPERATE THE MASSAGE CHAIR.....	8
I	CARE AND MAINTENANCEGUIDELINES.....	9
I	SPECIFICATIONS.....	10

IMPORTANT SAFETY INSTRUCTIONS

OPERATING ENVIRONMENT

- Do not use or store the massage chair under damp conditions:
- i.e. bathroom, sauna, spa, to avoid any mold development, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors
- Don't use the massage chair under dusty and temperature varying conditions.
- Don't place the massage chair near heater(s), stove or direct sunlight to avoid inner components overheating which may lead to electrical fires and shock.
- Make sure to place the massage chair on flat level ground surface
- The chair should be placed on carpet, wood or hard surface flooring
- Make sure the power plug is properly grounded to avoid shock and injury.
- Don't use the machine while smoking



FOLLOWING CONDITIONS REQUIRING CONSULT FROM A MEDICAL DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Bone disease
- While pregnant or menstruating
- On rest by the doctor
- Spinal condition including injuries associated to the back and spine
- Hypertension and cramping
- Arthritic bone disease
- Any implanted medical device such as pacemaker
- Suffering from internal conditions (i.e. gastritis, enteritis, diarrhea, Zepatitis)
- Serious blood circulation trouble
- Consciousness issues such as vertigo and other equilibrium matters

WARNING

Children under the age of 18 years old or with physical or mental condition are not permitted to use this product unattended by an able adult and with the permission of a medical doctor. These precautions are taken to prevent injury and damage to the product

IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

⊘ CAUTIONS

Non-authorized technicians are not permitted to dismantle or repair this product to avoid fire, electrical shock or serious injury



Do not clean or wipe down this product with anything wet or moist when this symbol is shown.



When storing or during maintenance, be sure to disconnect the power plug from socket



Do not detach the plug with wet hands

Make sure to turn the power switch to the "off" position and detach the power plug when cleaning or storing the product



During cleaning or when moving the chair, make sure to detach power cord from the wall and chair to avoid fire or electrical shock

Do not place the machine on a unstable, uneven and unlevelled surface to avoid damage, insulation failure or electrical shock



When storing the machine, make sure to disconnect power plug from the outlet and chair to avoid fire caused by electrical surges.

Clean power plug periodically to avoid buildup leading possible fire



Do not use it at any unsafe area or for any other purposes than personal use. Not for commercial use.



Do not store the chair in areas that are not climate controlled. Must be in room temperature with very low humidity.



Keep the power cord unobstructed at all time, i.e. not twisted, knotted, over bent or pulled and do not place anything heavy on it



IMPORTANT SAFETY INSTRUCTIONS

⊘ CAUTIONS

**Do not make any modifications or change components by an unauthorized person
Do not dismantle the machine by unauthorized person to avoid electrical shock,
short-circuit or fire**

**Do not store it near a fire, heater or flammable substances
It is one of the reasons of causing product damage.**



**Should any abnormality occur(i.e. burning smell), stop using immediately
and detach power plug from outlet.
please contact your local seller or authorized service agent for maintenance
to avoid electrical fire and injury.**



Only proper rated voltage must be used



**Stop using immediately if power cord or power plug is damaged or
not tightly fitted into the outlet or the base of the chair**



Be sure to keep children and pets away while operating to avoid any injury



Do not place the chair in direct sunlight at all times




**If the supply cord is damaged,it must be replaced by the manufacturer or an
authorized technician**



PARTS NAME

NAME OF EACH PART

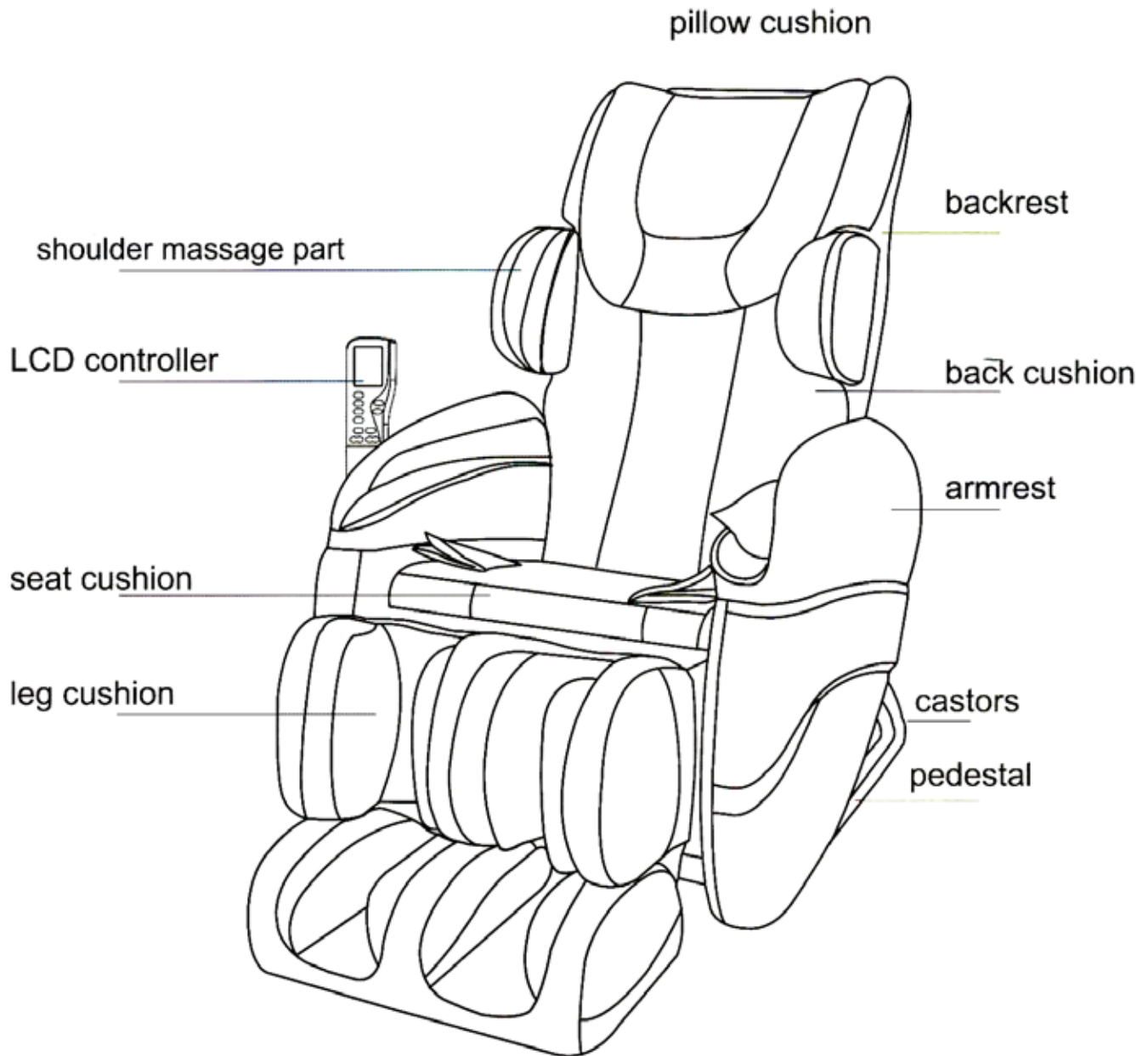


WARNING

1. Make sure the power cable and any air hoses are NOT twisted or bent. Do not place any object on top of the power cable.
2. Before plugging in the power cable into the outlet be sure the chair power switch on the back of the chair is in the off position.

INSTRUCTIONS ON CONTROL PANEL


HOW TO POWER ON AND POWER OFF

1. Connect the power line and switch the power button to the "ON" position

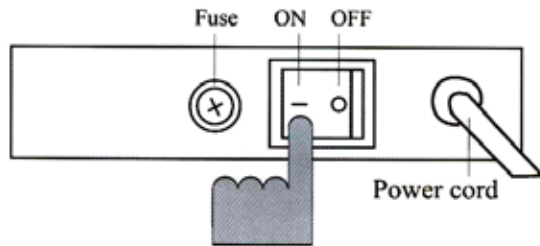

2 How to start the massage chair

When you are comfortably sitting in the chair, press the power button. Then select the massage program you would like to run. In the auto program it will run an initial body scan. The LCD remote will indicate what the chair is doing and the location of the rollers. Press the ON/OFF and it will close all programs, then the backrest and leg rest will return to default position.

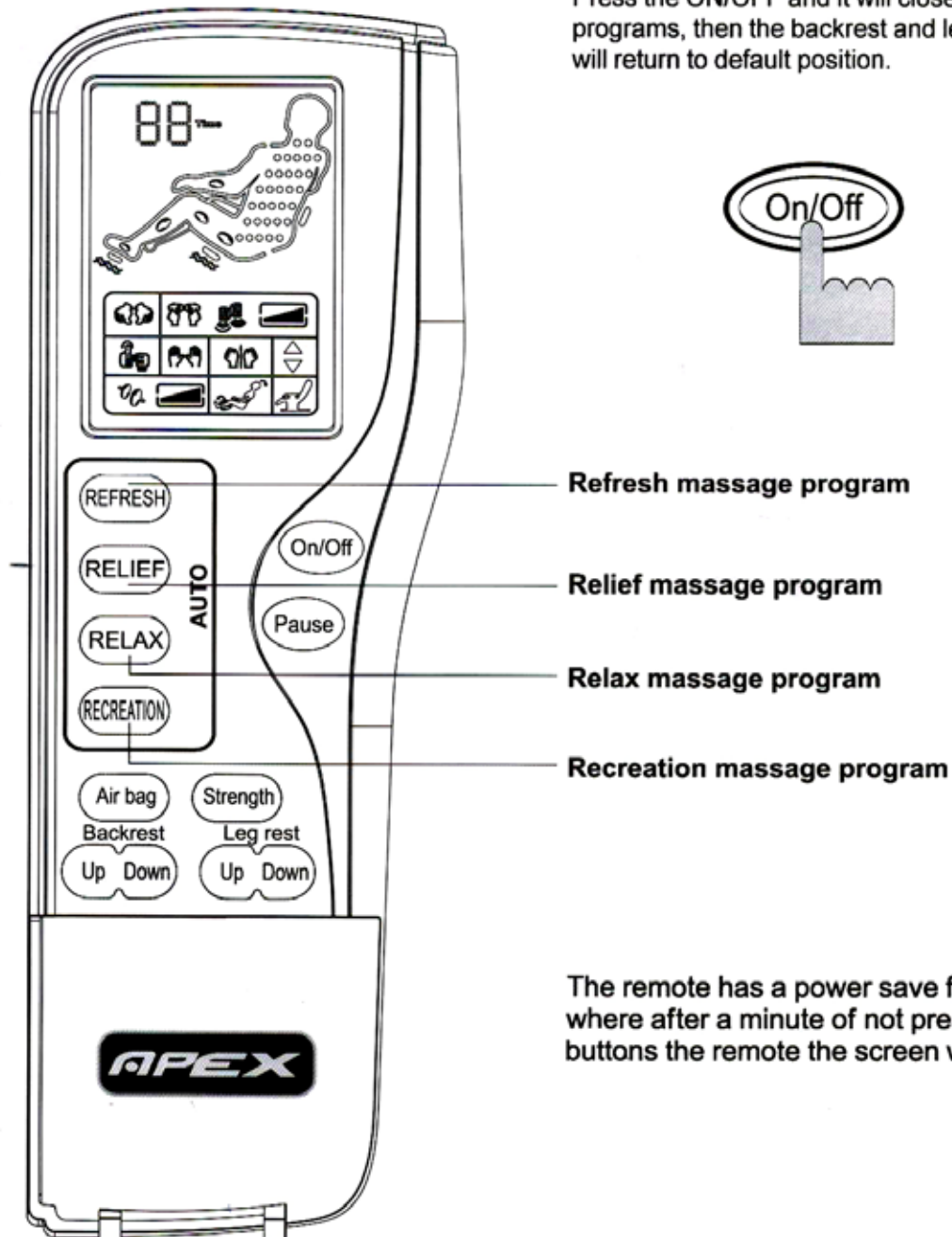

Refresh massage program


Relief massage program

Relax massage program

Recreation massage program

The remote has a power save feature where after a minute of not pressing any buttons the remote the screen will go dark


INSTRUCTIONS ON CONTROL PANEL

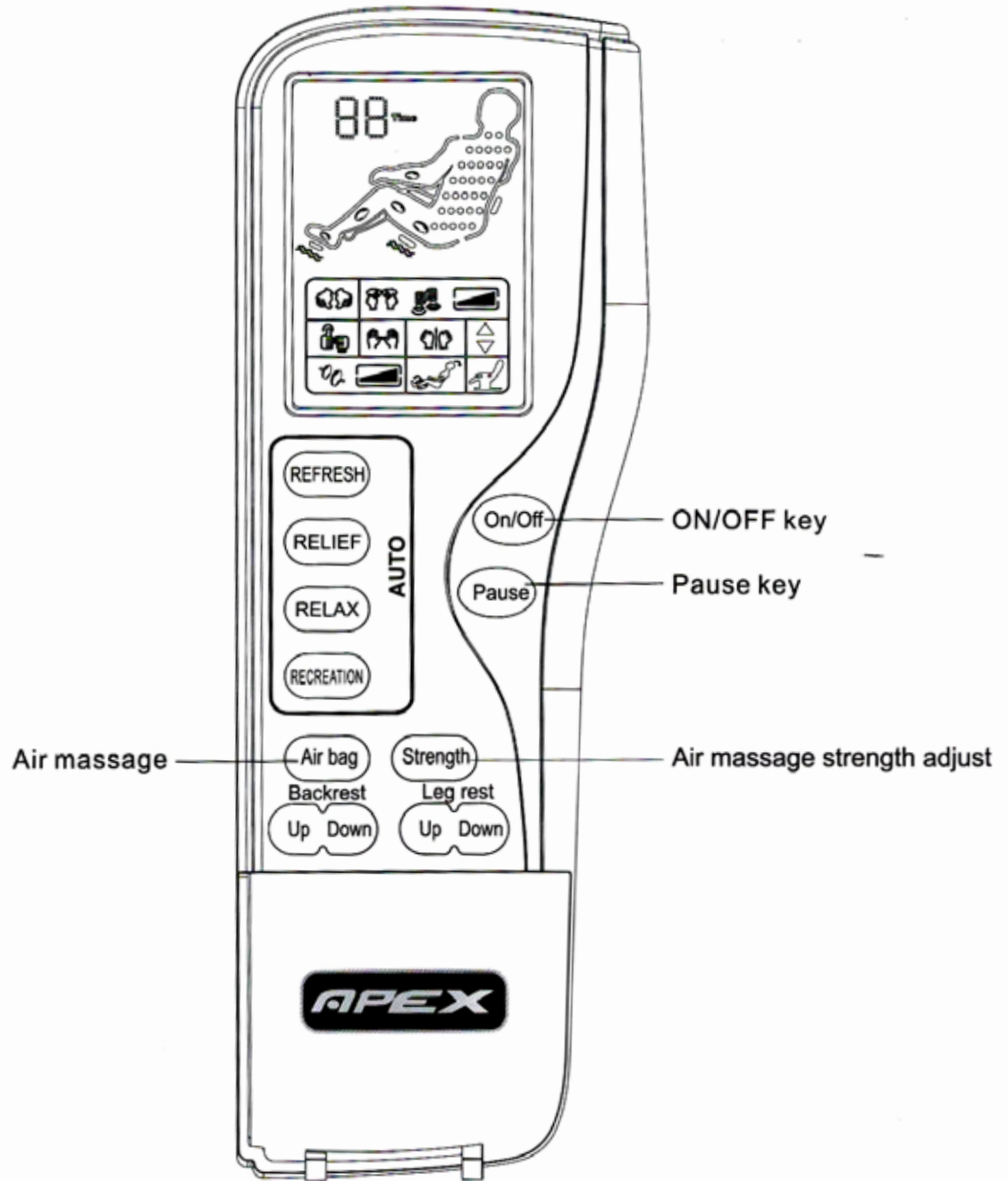


3. Adjust the backrest and legrest

By pressing the backrest and legrest key it will adjust the angle, reclining or bring you back to an upright position.

Backrest

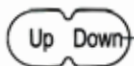


Legrest

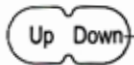

Press this key backrest up/down and legrest up/down are adjust the angle together

INSTRUCTIONS ON CONTROL PANEL

INSTRUCTIONS OF BUTTONS

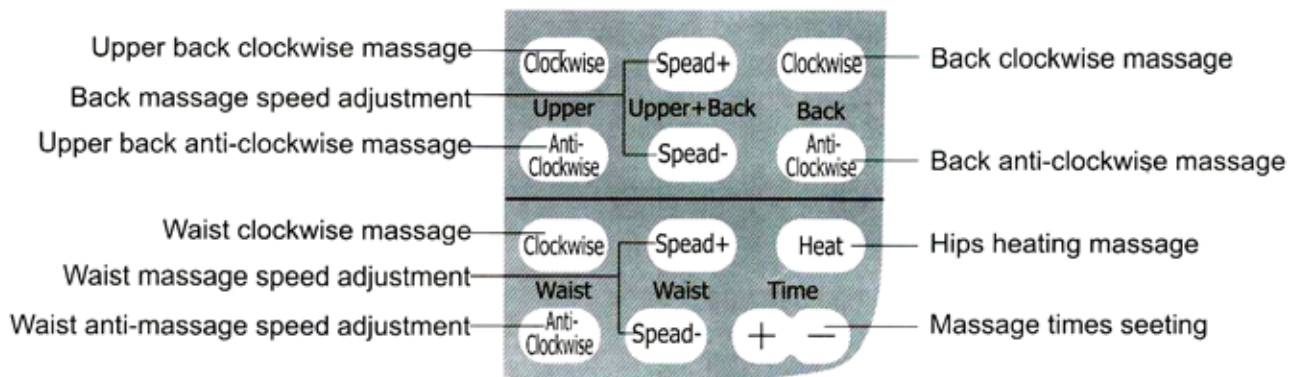

4. Difference from stop and finishing a massage

End of a preset program: At the end of a massage program the chair will stay in the reclined position and will not move to its default setting.

Forced stop: When you stop the massage in the middle of the program it will stay in the position it was stopped at. If you press the stop button one more time it will recline up to its default position.

Avoid the legrest damage, please lie the legrest to the initial condition when you stop massage.

5. Start the massage

There are 4 Auto massage programs

REFRESH, RELIEF, RELAX AND STRETCH

Auto massage program has different massage technique:

Thumb—press, Kneading, Tapping and Air press massage.

KEYS	DESCRIPTION	MAIN FUNCTION
REFRESH	Refresh massage program	kneading massage, air massage and tapping massage. You can adjust the angle during massage.
RELIEF	Relief massage program	Back massage and leg massage. You can adjust the angle during massage.
RELAX	Relax massage program	Slow massage and air massage. You can adjust the angle during massage.
RECREATION	Recreation massage program	Back and leg stretch massage. The angle can NOT be adjusted during massage.








All of the auto massage program, the time frame is set for 15 minutes, air bag strength is set to a default medium intensity. also you can press the STRENGTH button to adjust the strength of air bag.

All the auto massage program is stand alone program.

HOW TO OPERATE

INSTRUCTIONS OF PARTS AUTO MASSAGE PROGRAM

CLOCKWISE AND ANTI-CLOCKWISE


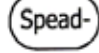


KEYS	DESCRIPTION	MAIN FUNCTION
 Upper Anti-Clockwis	Upper back clockwise massage	Massage speed & directional adjustable
	Upper back counter-clockwise massage	Massage speed & directional adjustable
 Upper+Back 	Upper back or whole back massage speed increase	3 speed adjustments for the massage
	Upper back or whole back massage speed decrease	3 speed adjustments for the massage
 Back Anti-Clockwis	Whole back clockwise massage	Massage speed adjustable
	Whole back anti-clockwise massage	Massage speed adjustable
 Waist Anti-Clockwis	Waist forward massage	Massage speed adjustable
	Waist reverse massage	Massage speed adjustable
 Waist 	Waist massage speed increase	3 levels of massage speed can be selected
	Waist massage speed decrease	3 levels of massage speed can be selected

The timing of auto massage program is 15 minutes, the default speed is set to medium, which can be adjusted by pressing the SPEED button.

HOW TO OPERATE

FUNCTION MESSAGE

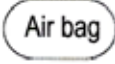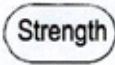AIR BAG AND TIMER

KEYS	DESCRIPTION	MAIN FUNCTION
	Air bag massage	This button turn on and off the air bag massage.
	The adjust key for air bag strength	Choose from 5 levels of air massage intensity.
	Time plus or minus key	This allows for the time of the massage to be reduced or increased.
	Lumbar Heat Massage	This button will activate the heating pad in the lumbar area

All the above is stand alone program that can be independently activated.

CARE AND MAINTENANCE

STORAGE

- Do not remove the backrest form the chair nor touch the internal mechanical parts of massage chair.
- Please follow these guidelines to avoid personal injury.
- Do not remove the internal parts of the chair nor repair the massage chair by yourself . If tampered with you are possibly subject to electrical shock and serious injury.
- Please pull out the power cord when you are cleaning the chair, if not it maybe subject to electric shock and serious injury. Do not use any wet clothes to wipe down the chair.
- The massage chair should be kept clean all the time , and set in the dry place at room temperature.
- If the massage chair not used for a long period of time ,please cover with a cloth to avoid dust getting into the chair. Avoid direct sunlight.

THE MAINTENANCE OF MASSAGE CHAIR

- When you clean the massage chair please use neutral detergent like vinegar water, and tightly ring out any liquids from the cloth.
- If the backrest cushion , pillow cushion and legrest cushion are dirty ,please use neutral detergent or vacuum out the dust and dirt.
- The plastic piping should be kept clean and free of any bending to allow it to operate properly.
- Do not use abrasive cleaning agents like benzene or thinner to clean the controller, this may result in malfunction and destruction of the remote.
- Do NOT use any aerosol cleaners on the chair. It may cause electrical fires. No Alcohol based agents should be used on the chair
- Please use dry cloth cleaning to avoid damage.

THE MAINTENANCE OF CONTROLLER

- Avoid any source of heat such as heaters, gas and ignitions.
- Keep dry at all times and free of dust.
- Please clean with a soft and dry cloth.
- Do NOT use and type of alcohol or benzene based cleaners which may ruin the remote

CARE AND MAINTENANCE

TROUBLESHOOTING

Please inspect the massage chair frequently for any tears, damaged wires. Particularly the backrest and behind the back pad. If there are any tears on the lining of the backrest or the back pad, stop using the chair and contact the retailer or supplier to have the problem resolved. This may cause serious injury or damage to the chair.

SITUATIONS	SOLUTION
The massage chair doesn't work, ON/OFF Restore key does not light	1. Check whether the plug is well connected and up and the power cord has no visible damage 2. Check whether the switch on the massage chair is turned to [I] on position back of the chair
Are you hearing any of the following mechanical sounds?	If you hear any of the following sounds, it is the normal function of the chair. 1. When Tapping and kneading a motorized and clicking sound will occur. If louder than normal please contact the manufacturer. 2. Massage wheels rubbing against the material. 3. Inflating and exhaust sound during air massage 4. Motorized sound when adjusting the recline 5. Slightly elevated motorized sound when user is above 200 lbs.

EMERGENCY TREATMENT

1. If the chair freezes unexpectedly or is not responding the remote, turn off the chair from the back power switch and turn back on to reset the chair.

2. If the power line is broken, please notify our customer service department for replacement part.
Warning: Only certified technicians are permitted to use the chair to avoid serious injury and damage to the chair. If tampered with the warranty we be voided.

SPECIFICATIONS

Description	Massage Chair
Model	AP-Prime II
Rated supply	110V~
Rated Frequency	60Hz
Rated Power	150W
Rated working hour	15 min