




AP-Odyssey

INSTRUCTION MANUAL

MASSAGE CHAIR

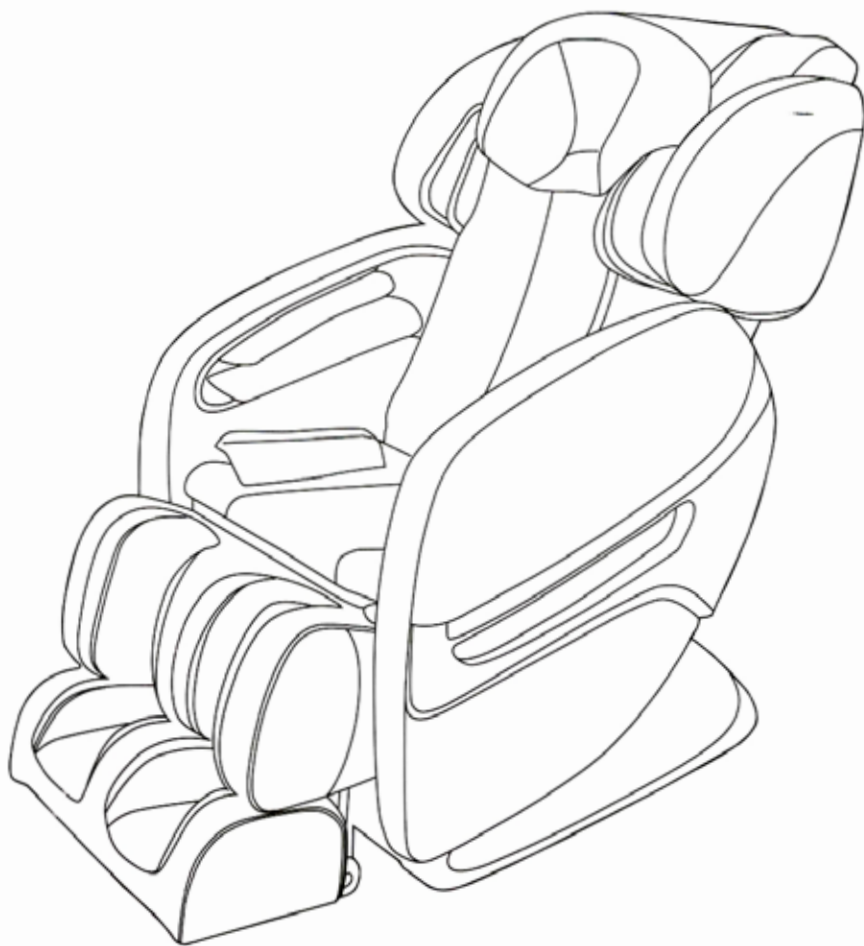

BE SURE TO READ THIS *INSTRUCTION MANUAL* CAREFULLY BEFORE OPERATING THE CHAIR. PLEASE KEEP THIS MANUAL IN SAFE KEEPING FOR FUTURE REFERENCE.

CONTENTS

Thanks for purchasing our massage chair, hopefully the product will provide will meet your massage needs. For service issues, please contact your local seller or manufacturer, never dismantle or repair this product by an unauthorized technician to avoid any injury or product damage.

I	IMPORTANT SAFETY INSTRUCTIONS	1-3
I	PARTS NAME.....	4
I	INSTRUCTIONS ON FUNCTIONS OF THE CONTROL PANEL	5-7
I	HOW TO OPERATE THE MASSAGE CHAIR.....	8
I	CARE AND MAINTENANCE GUIDELINES.....	9
I	SPECIFICATIONS.....	10

IMPORTANT SAFETY INSTRUCTIONS

OPERATING ENVIRONMENT

- Do not use or store the massage chair under damp conditions: i.e. bathroom, sauna, spa, to avoid any mold development, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors
- Don't use the massage chair under dusty and temperature varying conditions.
- Don't place the massage chair near heater(s), stove or direct sunlight to avoid inner components overheating which may lead electrical fires and shock.
- Make sure to place the massage chair on flat level ground surface
- The chair should be placed on carpet, wood or hard surface flooring
- Make sure the power plug is properly grounded to avoid shock and injury.
- Don't use the machine while smoking



FOLLOWING CONDITIONS REQUIRING CONSULT FROM A MEDICAL DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Bone disease
- While pregnant or menstruating
- On rest by the doctor
- Spinal condition including injuries associated to the back and spine
- Hypertension and cramping
- Arthritic bone disease
- Any implanted medical device such as pacemaker
- Suffering from internal conditions (i.e. gastritis, enteritis, diarrhea, hepatitis)
- Serious blood circulation trouble
- Consciousness issues such as vertigo and other equilibrium matters

WARNING

Children under the age of 18 years old or with physical or mental condition are not permitted to use this product unattended by an able adult and with the permission of a medical doctor. These precautions are taken to prevent injury and damage to the product.

IMPORTANT SAFETY INSTRUCTIONS

For safety reasons we have clearly indicate potential dangers if the chair is not used properly. Instructions have been marked in the following manners and should be strictly observed and followed.

⊘ CAUTIONS

Non-authorized technicians are not permitted to dismantle or repair this product to avoid fire, electrical shock or serious injury



Do not clean or wipe down this product with anything wet or moist when this symbol is shown.



When storing or during maintenance, be sure to disconnect the power plug from socket

Do not detach the plug with wet hands



Make sure to turn the power switch to the "off" position and detach the power plug when cleaning or storing the product

During cleaning or when moving the chair, make sure to detach power cord from the wall and chair to avoid fire or electrical shock



Do not place the machine on a unstable, uneven and unlevelled surface to avoid damage, insulation failure or electrical shock

When storing the machine, make sure to disconnect power plug from the outlet and chair to avoid fire caused by electrical surges.



Clean power plug periodically to avoid buildup leading possible fire



Do not use it at any unsafe area or for any other purposes than personal use. Not for commercial use.



Do not store the chair in areas that are not climate controlled.

Must be in room temperature with very low humidity.



Keep the power cord unobstructed at all time, i.e. not twisted, knotted, over bent or pulled and do not place anything heavy on it



IMPORTANT SAFETY INSTRUCTIONS

⊘ CAUTIONS

Do not make any modifications or change components by an unauthorized person

Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire



Do not store it near a fire, heater or flammable substances

It is one of the reasons of causing product damage.



Should any abnormality occur(i.e. burning smell), stop using immediately and detach power plug from outlet.

please contact your local seller or authorized service agent for maintenance to avoid electrical fire and injury.



Only proper rated voltage must be used



Stop using immediately if power cord or power plug is damaged or not tightly fitted into the outlet or the base of the chair



Be sure to keep children and pets away while operating to avoid any injury



Do not place the chair in direct sunlight at all times



If the supply cord is damaged,it must be replaced by the manufacturer or an authorized technician




WARNING

- 1.Make sure the power cord is not twisted, knotted nor anything heavy on top of it.
- 2.The ground wire must be connected properly
- 3.Before attaching the plug, make sure all switches are in the off position.
- 4.The massage mechanism will reset to the middle of the backrest which is the default setting.


PARTS NAME

NAME OF EACH PART

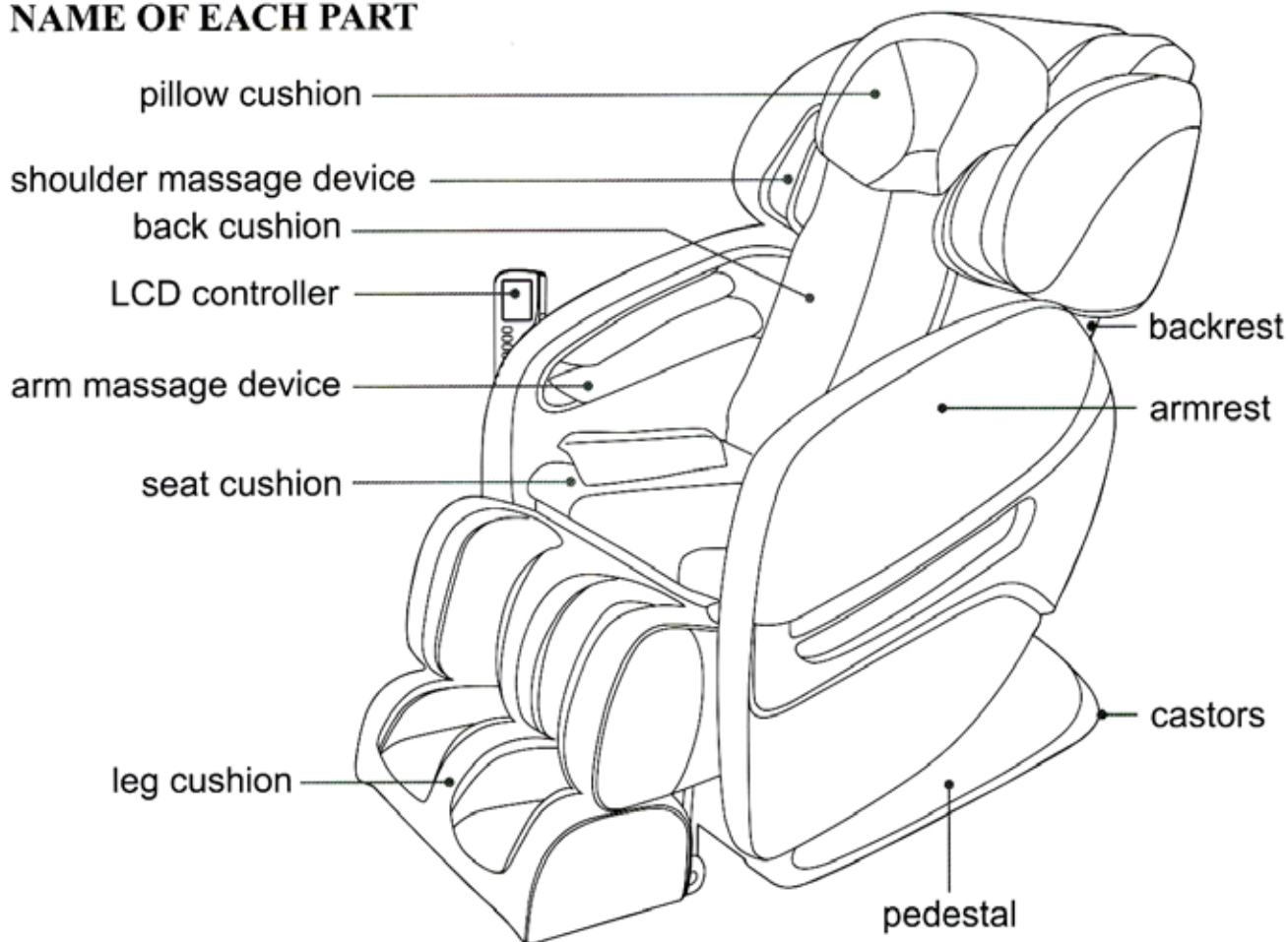

SHOULDER AND CAPE INSTALLATION

- ▶ First make sure the air pipe from mounting plate hole is pulled out to the full length, then insert the two bolts from the shoulder device into the backrest hole, tighten the screw, then connect the air pipes together and zipper up the cover.
- ▶ Make sure the air pipe at the top of backrest is fully extended, then connect the air pipe to the pipe from the pillow. Buckle the pillow to the backrest and zipper to the backrest.

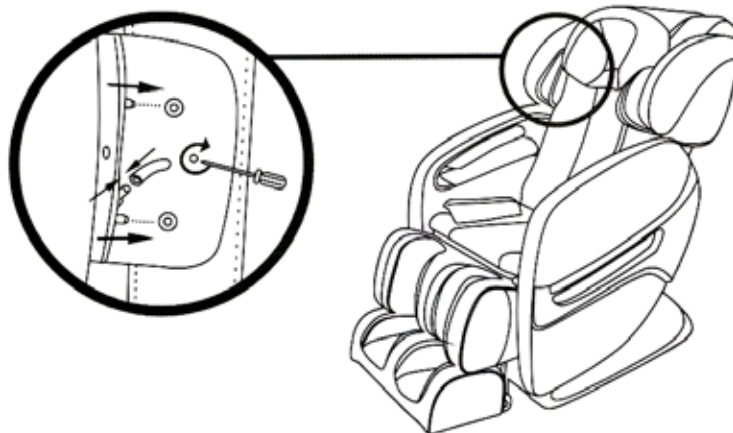


INSTRUCTIONS ON CONTROL PANEL

HOW TO POWER ON AND POWER OFF

1. Connect the power line and open switch



After connecting the power line, press the switch to [I] on position located at the back of the chair. The switch will light up indicating it has powered on.

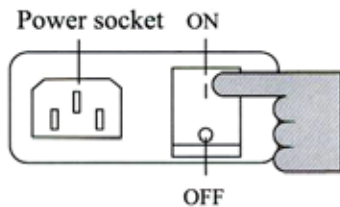

2 How to start the massage chair


When you are comfortably sitting in the chair, press the power button. Then select the massage program you would like to run. In the auto program it will run an initial body scan. The LCD remote will indicate what the chair is doing and the location of the rollers.

Press the ON/OFF and it will close all programs, then the backrest and leg rest will return to default position.

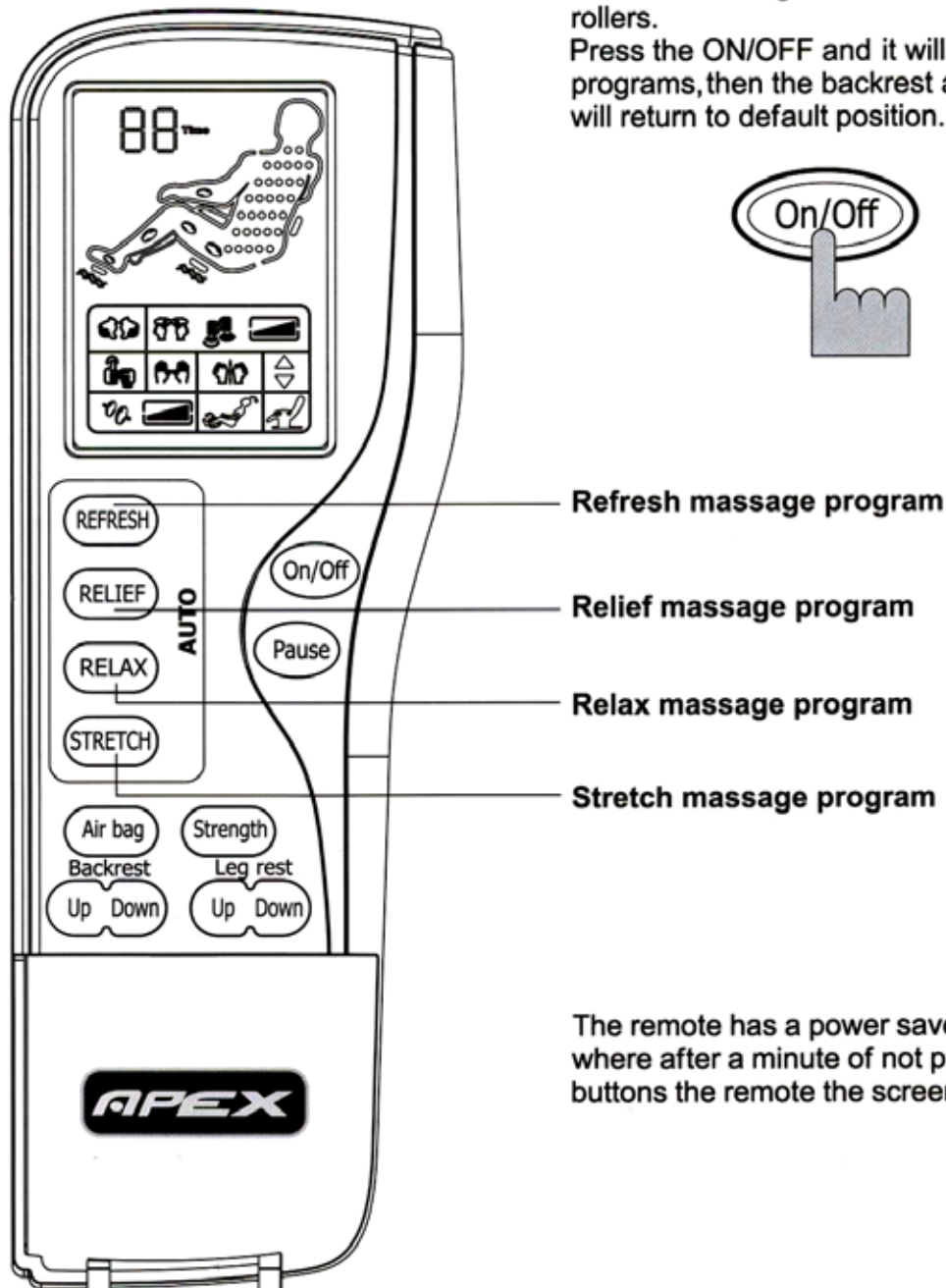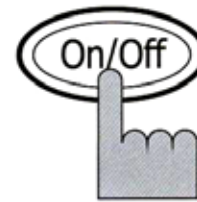

The remote has a power save feature where after a minute of not pressing any buttons the remote the screen will go dark.


INSTRUCTIONS ON CONTROL PANEL

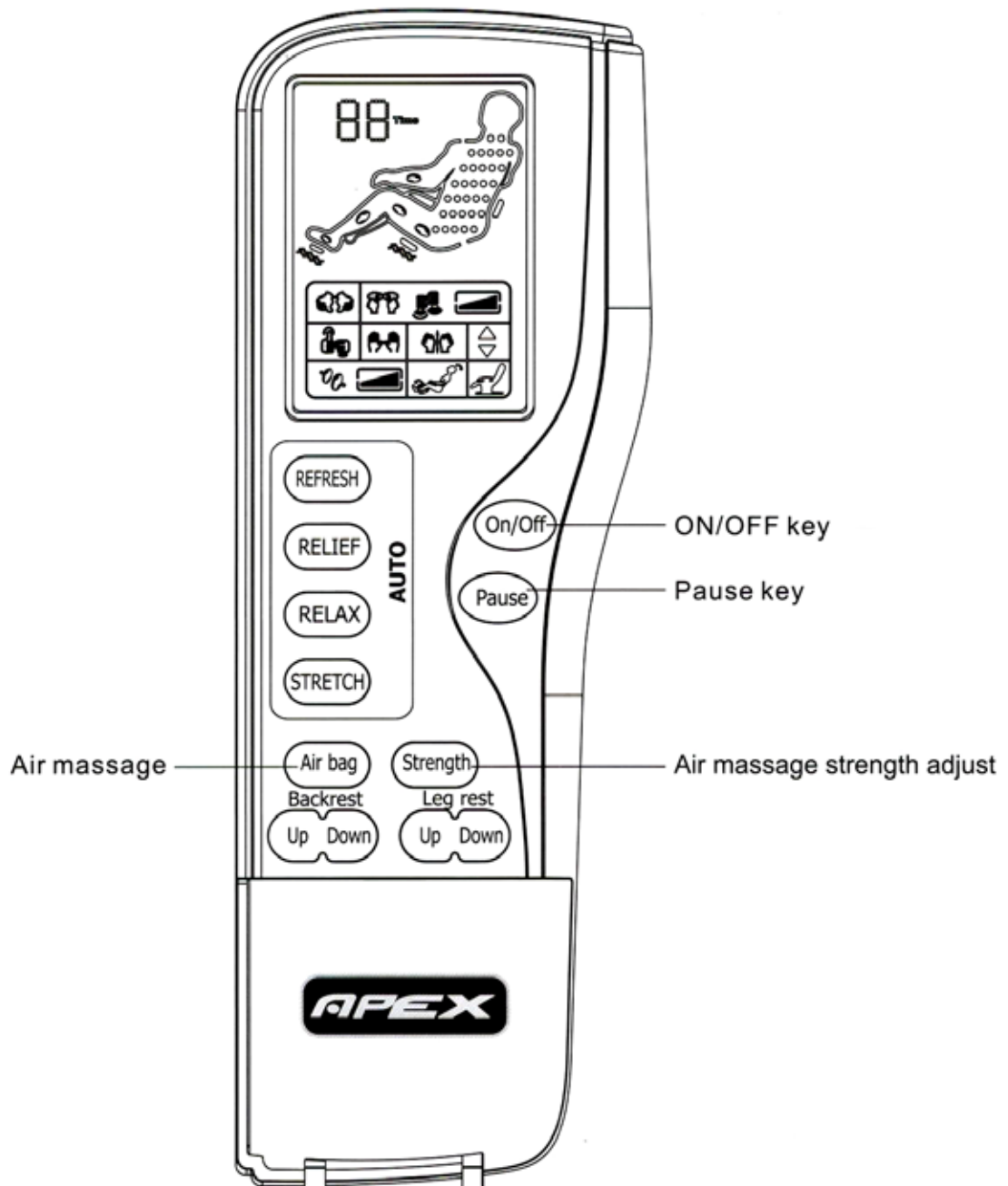

3. Adjust the backrest and legrest

By pressing the backrest and legrest key it will adjust the angle, reclining or bring you back to an upright position.

Backrest

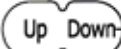
 Press this key backrest up/down and legrest up/down you will be able to adjust the backrest angle along with the footrest.

Legrest

 This key is only to control legrest and does not adjust the backrest angle.

INSTRUCTIONS ON CONTROL PANEL

INSTRUCTIONS OF BUTTONS

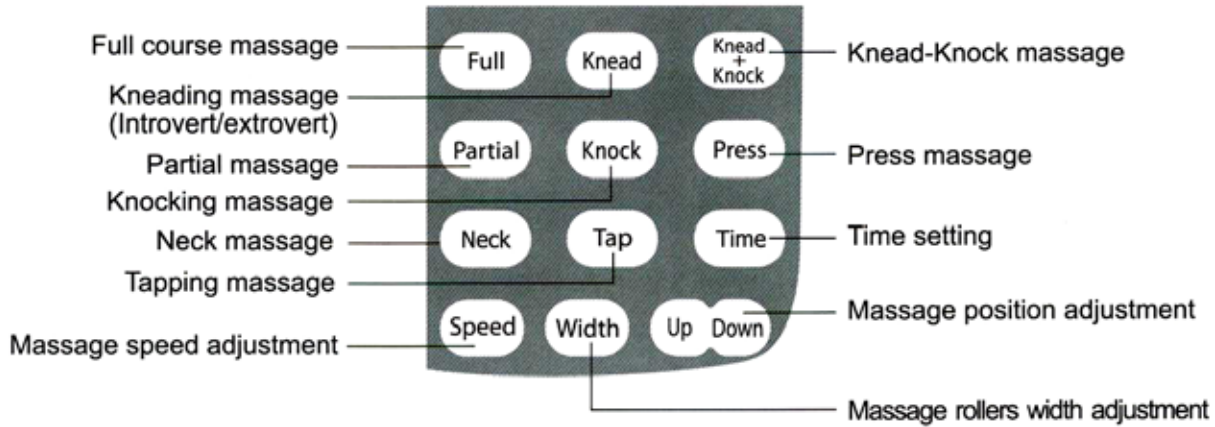

4. Difference from stop and finishing a massage

End of a preset program : At the end of a massage program the chair will stay in the reclined position and will not move to its default setting.

Forced stop: When you stop the massage in the middle of the program it will stay in the position it was stopped at. If you press the stop button one more time it will recline up to its default position.

Avoid the legrest damage, please lie the legrest to the initial condition when you stop massage.

5. Start the massage

There are 4 Auto massage programs

REFRESH, RELIEF, RELAX AND STRETCH

Auto massage program has different massage technique:

Thumb-press, Kneading, Tapping and Air press massage.

KEYS	DESCRIPTION	MAIN FUNCTION
REFRESH	Refresh massage program	kneading message, air message and tapping message. You can adjust the angle during message.
RELIEF	Relief massage program	Back message and leg message. You can adjust the angle during message.
RELAX	Relax massage program	Slow message and air message. You can adjust the angle during message.
STRETCH	Stretch massage program	Back and leg stretch message. The angle can NOT be adjusted during message.






All of the auto massage program, the time frame is set for 15 minutes, air bag strength is set to a default medium intensity, also you can press the STRENGTH button to adjust the strength of air bag.

All the auto massage program is stand alone program.

HOW TO OPERATE






INSTRUCTIONS OF USING THE MANUEL SETTING:

Air press massage, Kneading, Tapping, Press, Kneading and Tapping, Shiatsu, Kneading and Scratch Parts auto massage program aim to localize massage in the desired areas of the back.

KEYS	DESCRIPTION	MAIN FUNCTION
	Kneading massage	Back kneading massage program.
	Knocking massage	Back knocking massage program. You can adjust the width and speed during massage.
	Press massage	Back press massage program. You can adjust the width during massage.
	Kneading and Tapping massage	Kneading and Tapping massage program.
	Tapping massage	Back tapping massage program. You can adjust the width and speed during massage.

Please note that the manual mode and auto programs can not be used simultaneously. The timing of auto massage program is 15 minutes, the default speed of the massage is set to medium. You can press SPEED to adjust speed of the massage. The massage can be adjusted to a fixed point, partial or full back.

The key is related to parts auto massage program.

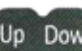KEYS	DESCRIPTION	MAIN FUNCTION
	The selection key of full back massage	Press it to start full back massage.
	The selection key of partial back massage	Press it to start partial back massage.
	The adjust key of massage speed	Choose five levels of speed.
	The adjust key of massage width	Choose three width: narrow, medium, width
	Change location of the massage heads	To adjust the position of massage heads, press UP or DOWN button to change its location.

HOW TO OPERATE

FUNCTION MESSAGE

AIR BAG AND VIBRATION

KEYS	DESCRIPTION	MAIN FUNCTION
Air bag	Air bag massage	To activate the air massage press the Air Bag button
Strength	The adjust key of air bag strength	Choose from five levels of air bag strength.
Neck	Neck massage	This will activate the massage that concentrates on the neck.
Time	Time setting	It's used to set massage duration.

CARE AND MAINTENANCE

STORAGE

- Do not remove the backrest from the chair nor touch the internal mechanical parts of massage chair. Please follow these guidelines to avoid personal injury.
- Do not remove the internal parts of the chair nor repair the massage chair by yourself. If tampered with you are possibly subject to electrical shock and serious injury.
- Please pull out the power cord when you are cleaning the chair, if not it maybe subject to electric shock and serious injury. Do not use any wet clothes to wipe down the chair.
- The massage chair should be kept clean all the time, and set in the dry place at room temperature.
- If the massage chair not used for a long period of time, please cover *with a* cloth to avoid dust getting into the chair. Avoid direct sunlight.

THE MAINTENANCE OF MASSAGE CHAIR

- When you clean the massage chair please use neutral detergent like vinegar water, and tightly ring out any liquids from the cloth.
- If the backrest cushion, pillow cushion and legrest cushion are dirty, please use neutral detergent or vacuum out the dust and dirt.
- The plastic piping should be kept clean and free of any bending to allow it to operate properly.
- Do not use abrasive cleaning agents like benzene or thinner to clean the controller, this may result in malfunction and destruction of the remote.
- Do NOT use any aerosol cleaners on the chair. It may cause electrical fires. No Alcohol based agents should be used on the chair
- Please use dry cloth cleaning to avoid damage.

THE MAINTENANCE OF CONTROLLER

- Avoid any source of heat such as heaters, gas and ignitions.
- Keep dry at all times and free of dust.
- Please clean with a soft and dry cloth.
- Do NOT use any type of alcohol or benzene based cleaners which may ruin the remote.

CARE AND MAINTENANCE

TROUBLESHOOTING

Please inspect the massage chair frequently for any tears, damaged wires. Particularly the backrest and behind the back pad. If there are any tears on the lining of the backrest or the back pad, stop using the chair and contact the retailer or supplier to have the problem resolved. This may cause serious injury or damage to the chair.

SITUATIONS	SOLUTION
The massage chair doesn't work, ON/OFF Restore key does not light up	1. Check whether the plug is well connected and the power cord has no visible damage 2. Check whether the switch on the massage chair is turned to [I] on position back of the chair
Are you hearing any of the following mechanical sounds?	If you hear any of the following sounds, it is the normal function of the chair. 1. When Tapping and kneading a motorized and clicking sound will occur. If louder than normal please contact the manufacturer. 2. Massage wheels rubbing against the material. 3. Inflating and exhaust sound during air massage 4. Motorized sound when adjusting the recline 5. Slightly elevated motorized sound when user is above 200 lbs.

EMERGENCY TREATMENT

1. If the chair freezes unexpectedly or is not responding the remote, turn off the chair from the back power switch and turn back on to reset the chair.
2. If the power line is broken, please notify our customer service department for replacement part.

Warning: Only certified technicians are permitted to use the chair to avoid serious injury and damage to the chair. If tampered with the warranty we be voided.

SPECIFICATIONS

Description	Massage Chair
Model	AP-Odyssey
Rated supply	110V~
Rated Frequency	60Hz
Rated Power	150W
Rated working hour	15 min