

Osaki

Item# : Healthmate-02

Shiatsu/Thermal Massage Bed

Instruction Manual




Preface

Thermal Health and massage Bed - the glad tidings to human health in the 21st century.

For the purpose of human health, with tenet of warmhearted service for customers and in the spirit of continual innovations, has developed the health and massage bed with multi-functions like health and massage. Thanks to its unique design, exquisiteness, novelty and the characteristics ahead of the times, it has been deeply favored by mass customers since its launch, fully embodying the health and massage effect of the bed.

Therefore, we decide to be oriented towards human and go up and up to make our efforts to go deeply into the health and Massage science and present our love to the health conditions of all mankind.

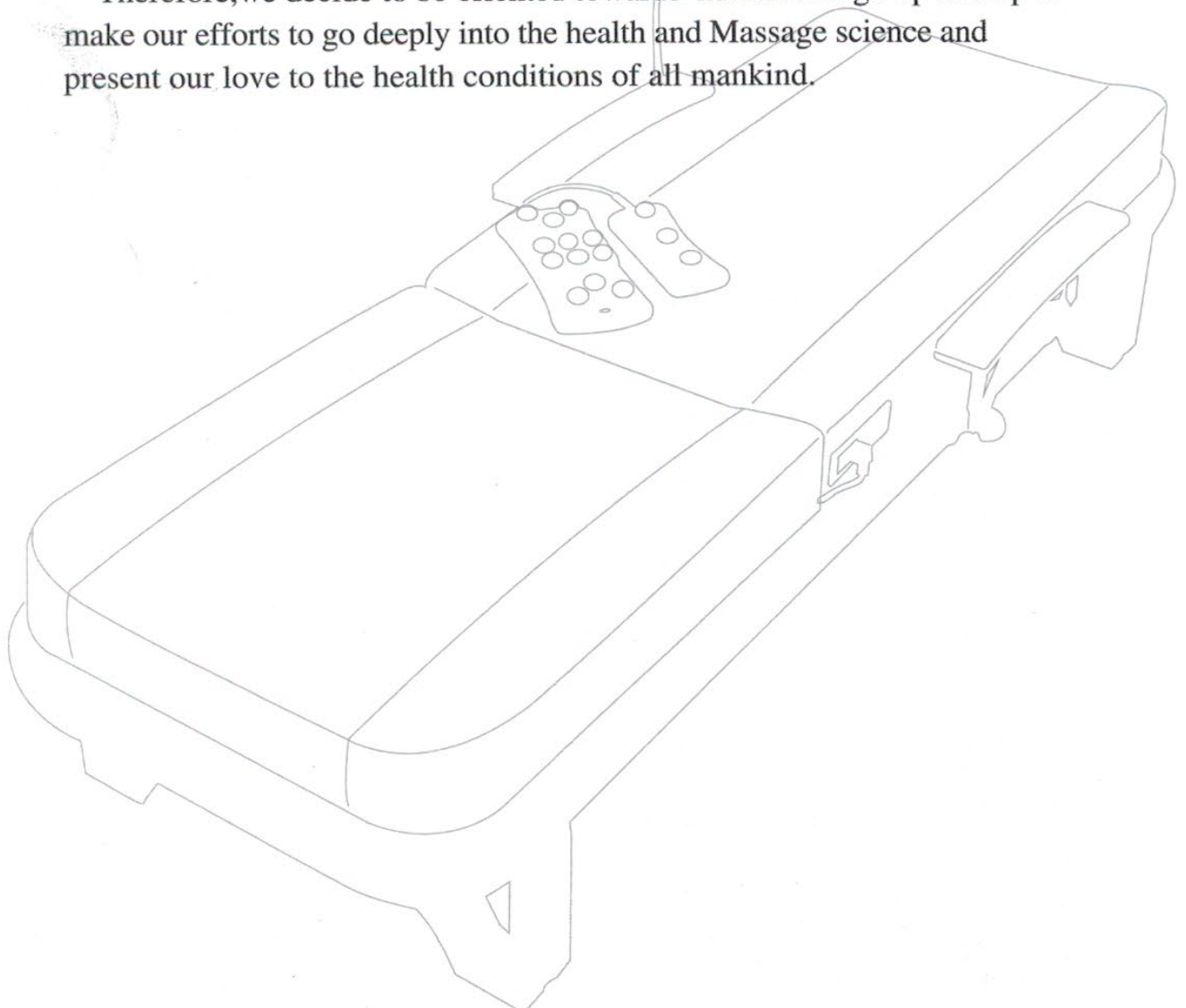

Contents

Contents



Characteristics -----	4
Structure, Performance and Working Principle	4
Technical Specifications -----	4
Names of the Parts of the Product -----	5
Diagram of the Use of the Remote Control --	6
Principle of the Product -----	7
Contraindications -----	7




Assembly Sketch -----	8-9
Product Installation -----	10
Trial and Precautions -----	10
Usage -----	10
Automatic mode to use -----	11
Various Modes -----	12-15



Safety Cautions -----	16
Notes for Use -----	16
Fee-based Services -----	17
A/S Service -----	18-19

Product Profile

IV. Names of the Parts of the Product

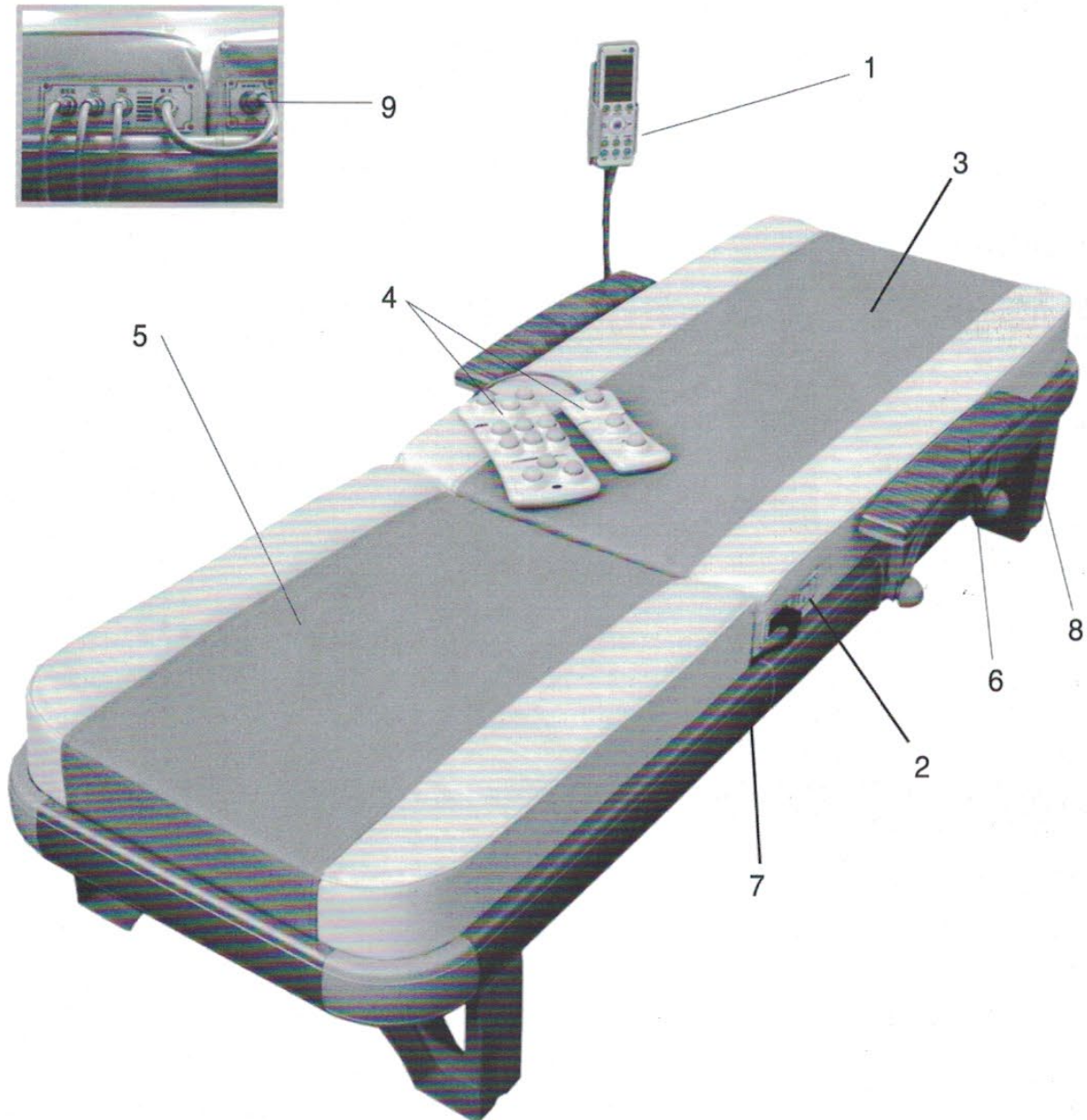


- | | | |
|--|---|----------------------|
| 1. Color-screen LCD remote control | 4. 3-ball and 12-ball external health heads | 7. Iron bed rack |
| 2. Connection between the remote control and external health heads | 5. Main part of the thermal health and Massage bed (back) | 8. Footboard |
| 3. Auxiliary part of the thermal health and Massage bed (legs) | 6. Armrests | 9. Power line socket |

Product Profile

Diagram of the Use of the Remote Control


Usage of the Remote Control

- | | |
|---------|--|
| Power | 1.Power Switch |
| AUTO | 2.The auto key has seven programs:Program1-Program7. Auto is set as default during start up, press the auto key and press ▼▲ key to select the mode and then press the "OK" key. |
| Stature | 3.Four modes of height. Press the "Stature" key to select 150cm/160cm/170cm/180cm. Press the ▼▲ key to select the height and then press the"OK" key. |
| Time | 4.Press the "Time" key to control massage time. Press the ▼▲ key to adjust the time. |
| | 5.Press the "I" Heating Element key and LCD displays "--"and press again to turn off. Press the "Pause"key first time: Both back roller and leg roller stop moving. Press the "Pause"key second time:back roller not move while leg roller start moving.Press the "Pause"key third time:back roller start moving while leg roller not move.Press the "Pause"key fourth time: Both back roller and leg roller start moving. |
| Manual | 6.Press the "Manual" key to select moving region. |
| Back | |
| Leg | |
| 12-ball | |
| 3-ball | |
- 7.Press the "Back"key and the LCD displays twinkle of the appointed part and then press the ▼▲ key to select the appropriate temperature. Finally press the "OK" key(press once for a rise or a drop of 5℃ in temperature, with the variation zone of 35℃-75℃)

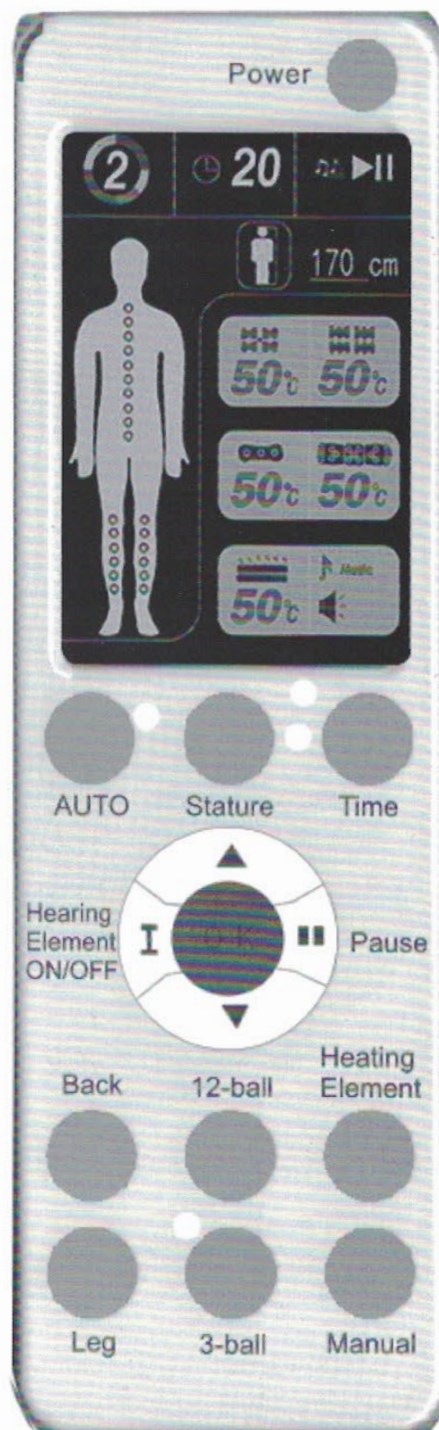

Product Profile

Description of LCD

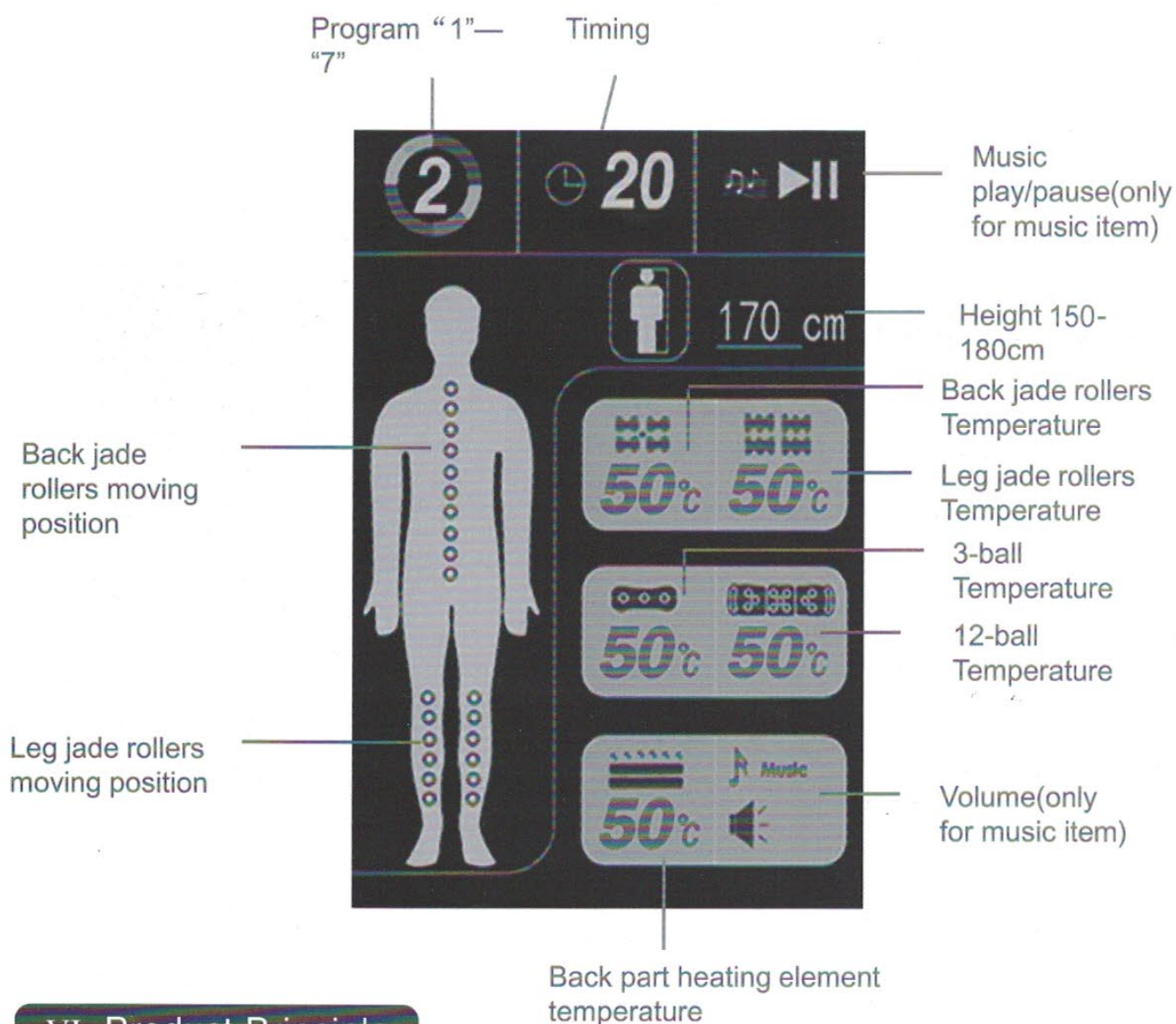

VI. Product Principle

Massage thermally to attain the effect of health and Massage.


VII. Contraindications

Pregnant and lying-in women, patients suffering from serious osteoporosis, malignant tumors, myocardial infarction, serious external wounds, high fever, serious hypertension and serious sequela of massive cerebral haemorrhage as well as extremely feeble persons are prohibited from using the product.




Product Installation

I. Please install as shown in the following figure.


Install bed frame

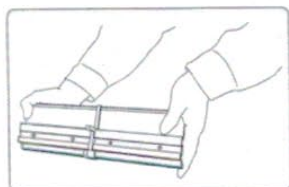



1、Put the plastic connection ring in the middle of aluminum internal connector.

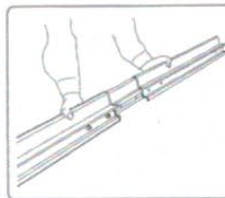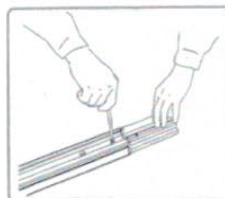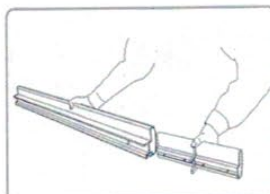

2、Equip one side of the aluminum frame with internal connector, then fix them with hexagon wrench and four screws (attention: the bigger screw hole on aluminum frame must be close to the plastic connection ring). Equip the other side of aluminum frame in the same way.





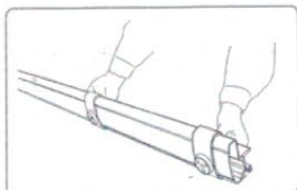


3、Put 2 pcs plastic armrest fixtures on each longer side of two aluminum frames

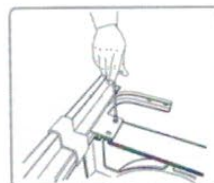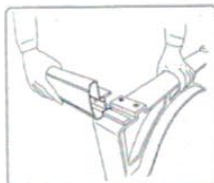

4、Insert two ends of one side aluminum frame into the bed legs, same way for the other side of aluminum frame.

5、Equip an iron bend to the four corners of bed frame respectively, and screw each holes on iron bend but do not tight up first.

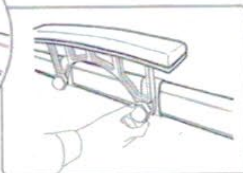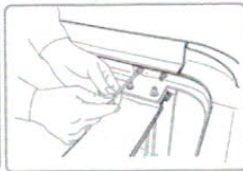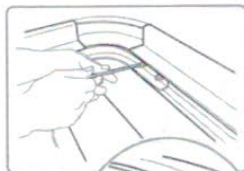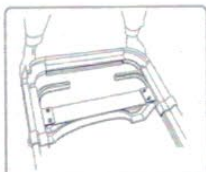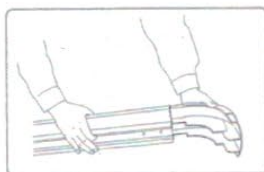


6、Insert the plastic bend to two sides of the shortest aluminum frame respectively, then insert two sides of aluminum frame with plastic bend to the long aluminum frame, at last tight up 6pcs screws on each iron bend to solid the bed frames and bed legs.(note: each set screw with 1pc M6 hexagon socket screw, 1pc spring washer and 1pc plain washer)

7、Each foot of armrest have three holes, fix the holes to plastic armrest fixtures with M4 size sunk screw, finished as the picture shows armrest. Each foot of armrest equip with plastic handwheels then fix the holder of remote controller and adjust the direction for user's convenient operation.


Product Installation

I. Please install as shown in the following figure.

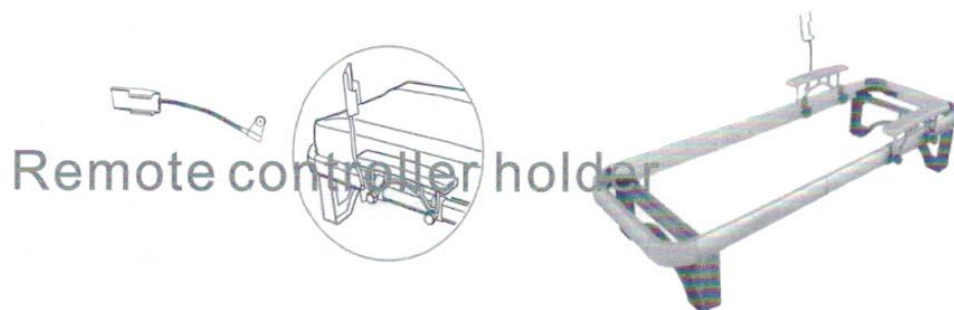


8. At last, check and tight up all screws on aluminum bed frame. Install bed frame work finished completely.

Install main part and leg part of jade massage bed

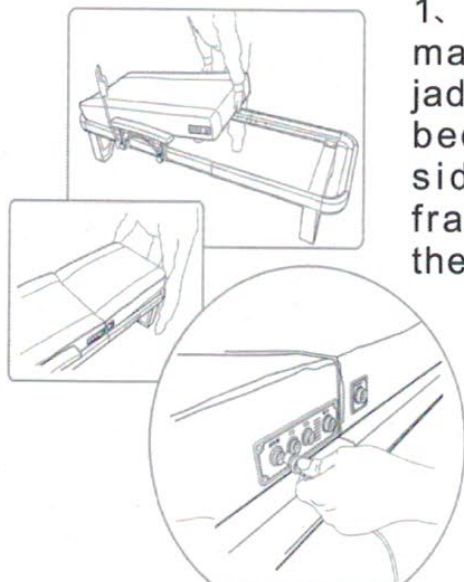


1. Place the main part of jade massage bed on long side of bed frame, place the leg part of

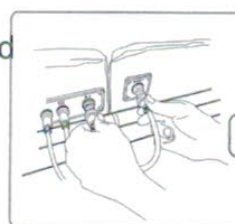

2. According the sticker indication on the side of jade massage bed, then connect the plug, remote controller, 3 balls and 12 balls jade handle equipment to the relevant sockets on bed (as marks on sticker of jade massage bed)

3. Whole of the work about install jade massage bed finished completely, please connect it with correct power

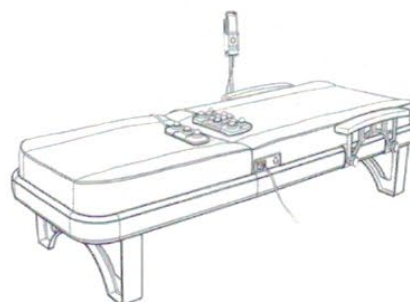


Product Installation & Usage

II. Product Installation

1. First place the main part and the auxiliary part on the bedstead and then connect the main part and the auxiliary part.
2. Get out of the power line, remote control, remote control bracket, 3-ball health head, 12-ball health head, power connecting lines of the main part and the auxiliary part and armrests.
3. After inserting the remote control into the remote control bracket, connect to the relevant sockets of the main part and fix the slots with nuts.
4. Connect the 3-ball health head and the 12-ball health head to relevant sockets of the main part in the same way and fix the slots with nuts.
5. Insert the other end of the power line into the slot of the thermal health and Massage bed and then connect it to 220V power supply.
6. If first connect to 220V power supply before plugging the thermal health and Massage bed in, it might cause a fault.

III. Trial and Precautions

1. Switch on power. (Note: switch off the main power supply in idle hours).
2. Press the power switch key at the top right corner of the remote control, press automatic key mode and then LCD displays  in a rolling way in the outer ring at the lower left. Press the "Confirmation" key, the internal health head starts to run. (Refer to the usage of the product for the specific usage.)
3. Be sure to connect the remote control before connecting power.
4. Temperature can be set at 45°C–55°C as users consider the most appropriate.
5. Prepare cotton clothes as thin as bedgown and unstrap and untie neckties.


IV. Usage

1. Lie on your back with your head towards the main part (back) of the thermal health and Massage bed and feet towards the auxiliary part (leg) of the bed.
2. In the automatic mode, massage of the abdomen or the hand with the 3-ball or 12-ball external health head can get a perfect effect.
3. In the supine posture with the face towards the main part of the thermal health and Massage bed, massage of the abdomen with the manual key (note: massage of the chest is not allowed) can also achieve a good effect.
4. One person uses in the basic way.
5. Two persons can use the product at same time. The other person puts the 3-ball external health head under the buttocks and massage all the parts of the body with the 12-ball external health head.
6. The "Automatic" or "Manual" mode can be used over and over to massage the parts needing to be massaged.



Various Modes

Automatic mode to use



The auto key has seven programs/modes: Program1- Program7. Auto is set as default during start up, press the auto key and press ▼▲ key to select the mode and then press the "OK" key. When turning on the power, the default mode is auto mode1.

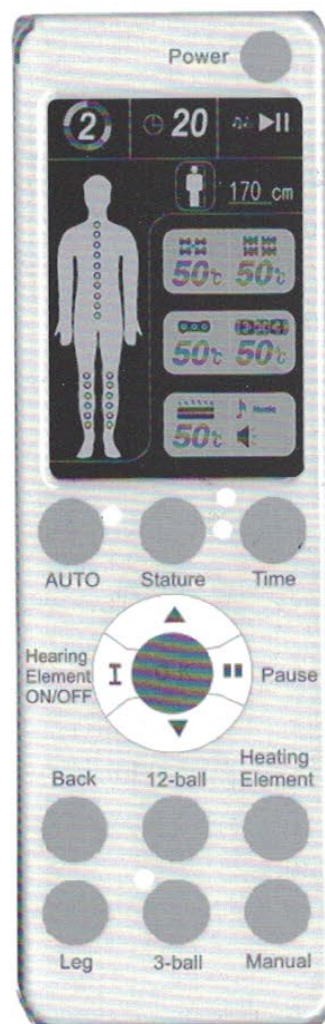


Various Modes

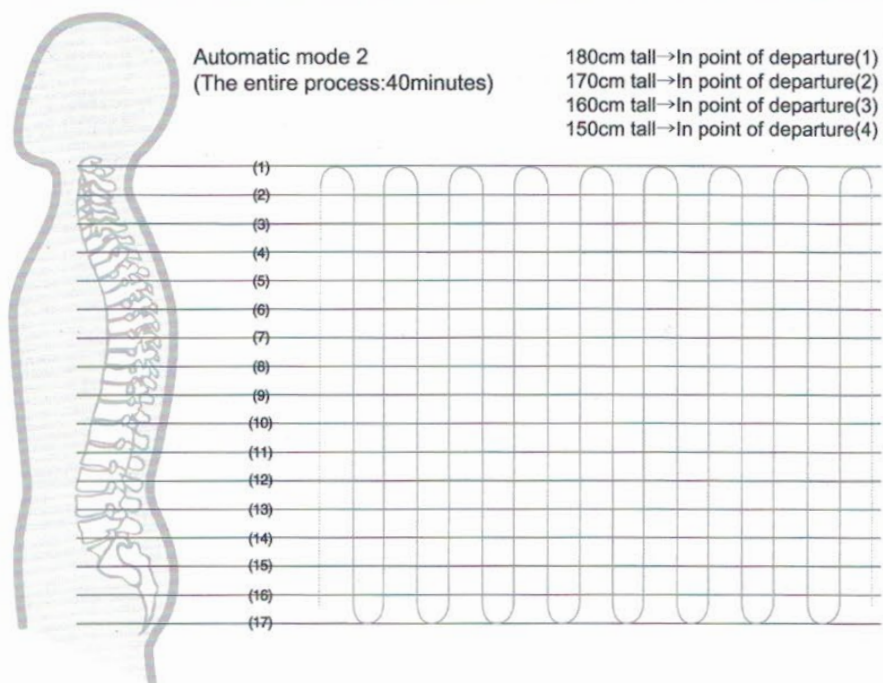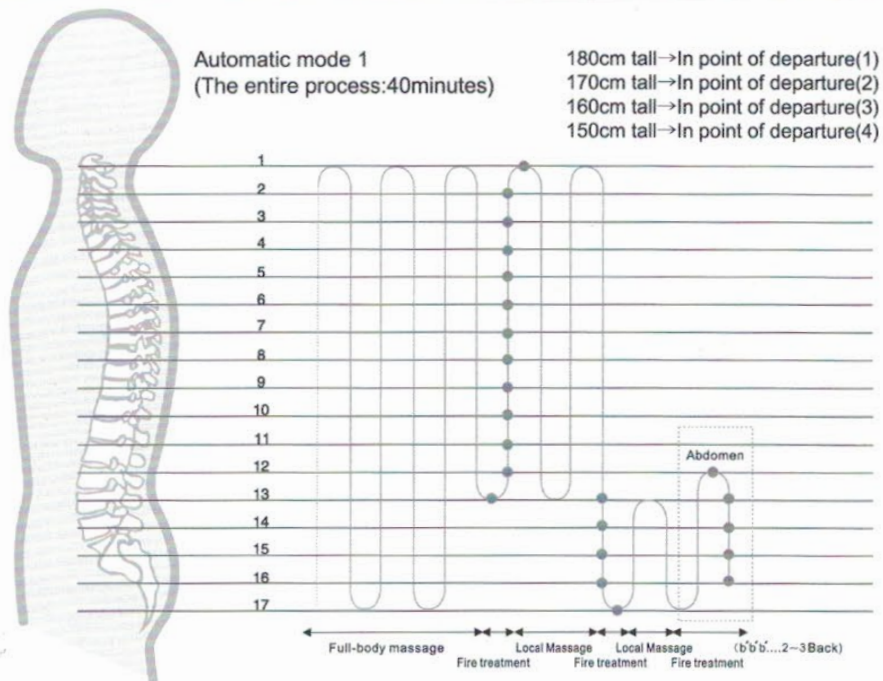

Various Modes

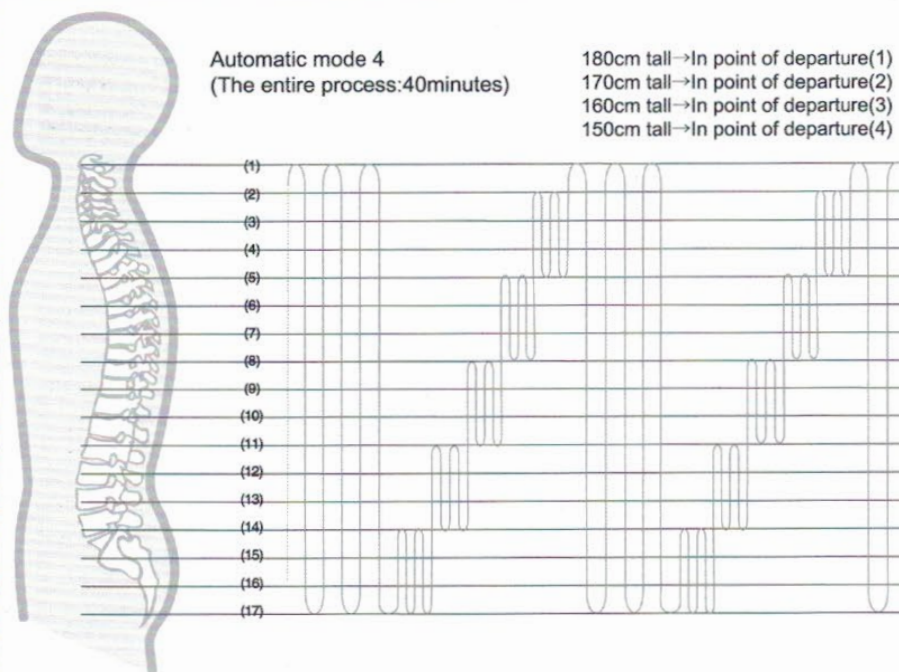
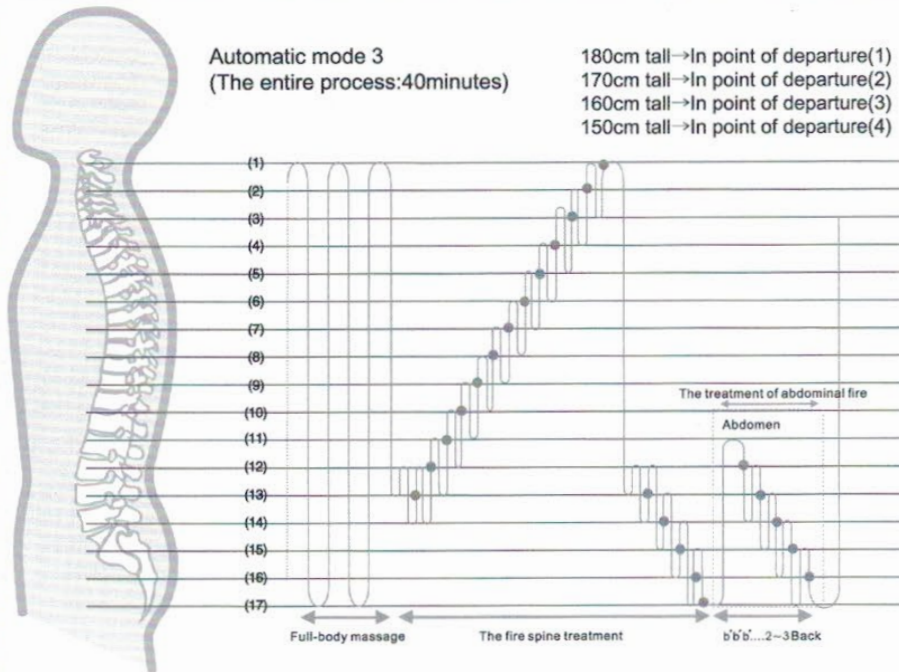


Various Modes

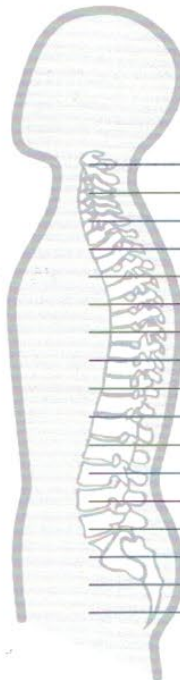


Automatic mode 5
(The entire process:40minutes)

180cm tall→In point of departure(1)
170cm tall→In point of departure(2)
160cm tall→In point of departure(3)
150cm tall→In point of departure(4)

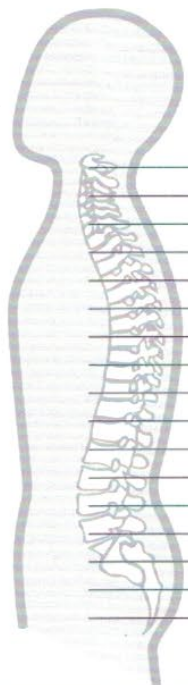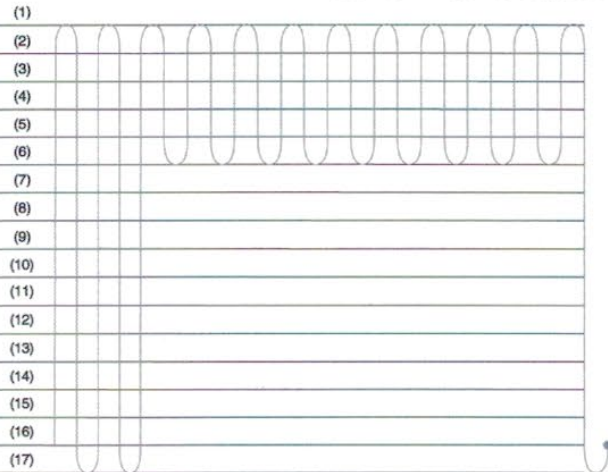

Automatic mode 6
(The entire process:40minutes)

180cm tall→In point of departure(1)
170cm tall→In point of departure(2)
160cm tall→In point of departure(3)
150cm tall→In point of departure(4)

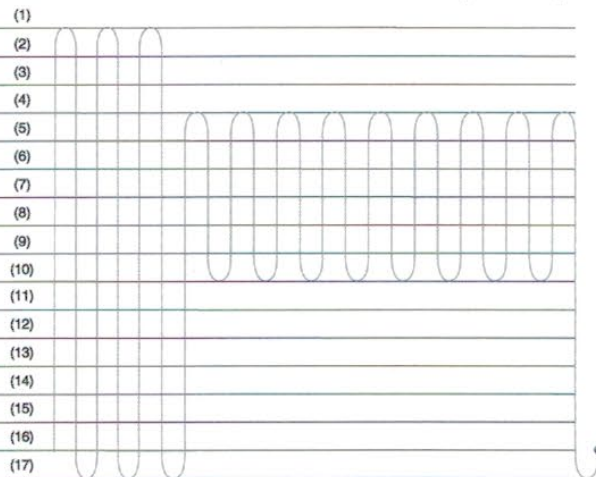


Various Modes



Automatic mode 7
(The entire process:40minutes)

180cm tall→In point of departure(1)
170cm tall→In point of departure(2)
160cm tall→In point of departure(3)
150cm tall→In point of departure(4)

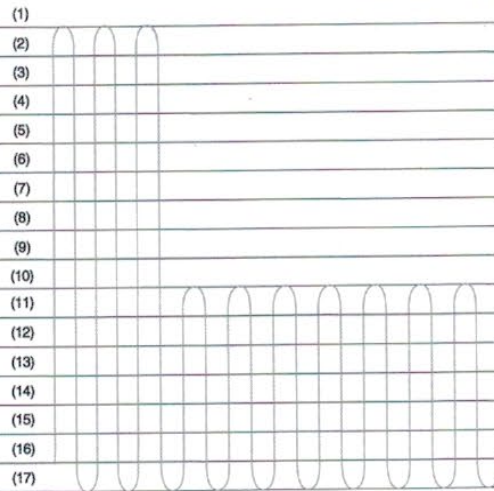

Precautions

I. Safety Cautions

1. When use comes to an end or the product is idle for a long time, you must switch off the power of the remote control and then unplug.
2. In case of any abnormality in the thermal health and Massage bed, please unplug and contact the after-sale service center.
3. Please don't bend the power plug or put any heavy article on it at random.
4. Please don't use any damaged power line plug or socket.
5. Please don't unplug violently or touch the power plug with your wet hand.
6. Please don't insert many plugs in one socket.
7. Please unplug in thundering and lightning.
8. Except the health head cover offered by our company, keep all the articles (fibre products, bedding, etc.) that may cause a fire away from the health heads.
9. As too high temperature may lead to scald, please set use temperature as (45°C-55°C) appropriate.
10. In installing the product, please make sure about 220V power supply.
11. In case of any trouble with the thermal health and Massage bed, you must contact the personnel of the after-sale service center to maintain or replace instead of disassembly and assembly by yourself.

II. Notes for Use

● Notes before Use

- Whether the power line and the socket are connected;
- Whether the socket screws of the main part and the health heads are fixed;
- Whether the power switch and the remote control are in normal operation;
- Whether temperature is displayed accurately;
- Whether the temperature of the health heads is within the appropriate use temperature range (45°C-55°C);
- Whether the health heads are abnormal in appearance;
- Whether the cover of the health heads is damaged.

● Notes after Use

- Check whether the remote control LCD is turned off
- Check whether the power switch is OFF
- Check whether the power line is disconnected from the host

● Notes in Use

- In using, please cover the health heads with the cover to avoid direct touch of skin (with the danger of scald).
- In order to avoid eyeball (pupilla) hurt, please don't put a health head on your eyeballs in using.
- Children or retardates should use the product accompanied by adults (in idle hours, please keep the product away from children). Users suffering from skin allergy should use at a lower temperature (to avoid the danger of scald).
- Please comply with the time of service and usage (use the product three-four times per day and please don't use any other article except for the articles provided by the manufacturer).

Fee-based Service

III. Fee-based Service

- Service is charged for the problem not arising from the product quality (caused by customers' improper keeping instead of the quality problem of the product itself).
- Service is charged for faults caused by customers' improper use.
- Service is charged when the consumable parts are completely abraded.
- In case of a fault caused by an event of force majeure, the fees arising therefrom should be charged.
- CureSense product is guaranteed for three years. When the product is used for business purpose and in the places of frequent service, the warranty period will be shortened to one year.


IV. Record of Events

[illegible]

* Note: This product has been strictly tested before delivery and our company reserves the right to modify the technical performance and device of the product. All the technical parameters are subject to change without further notice and the changed parameters should prevail.

A/S Service

I. Massage bed electric function diagram

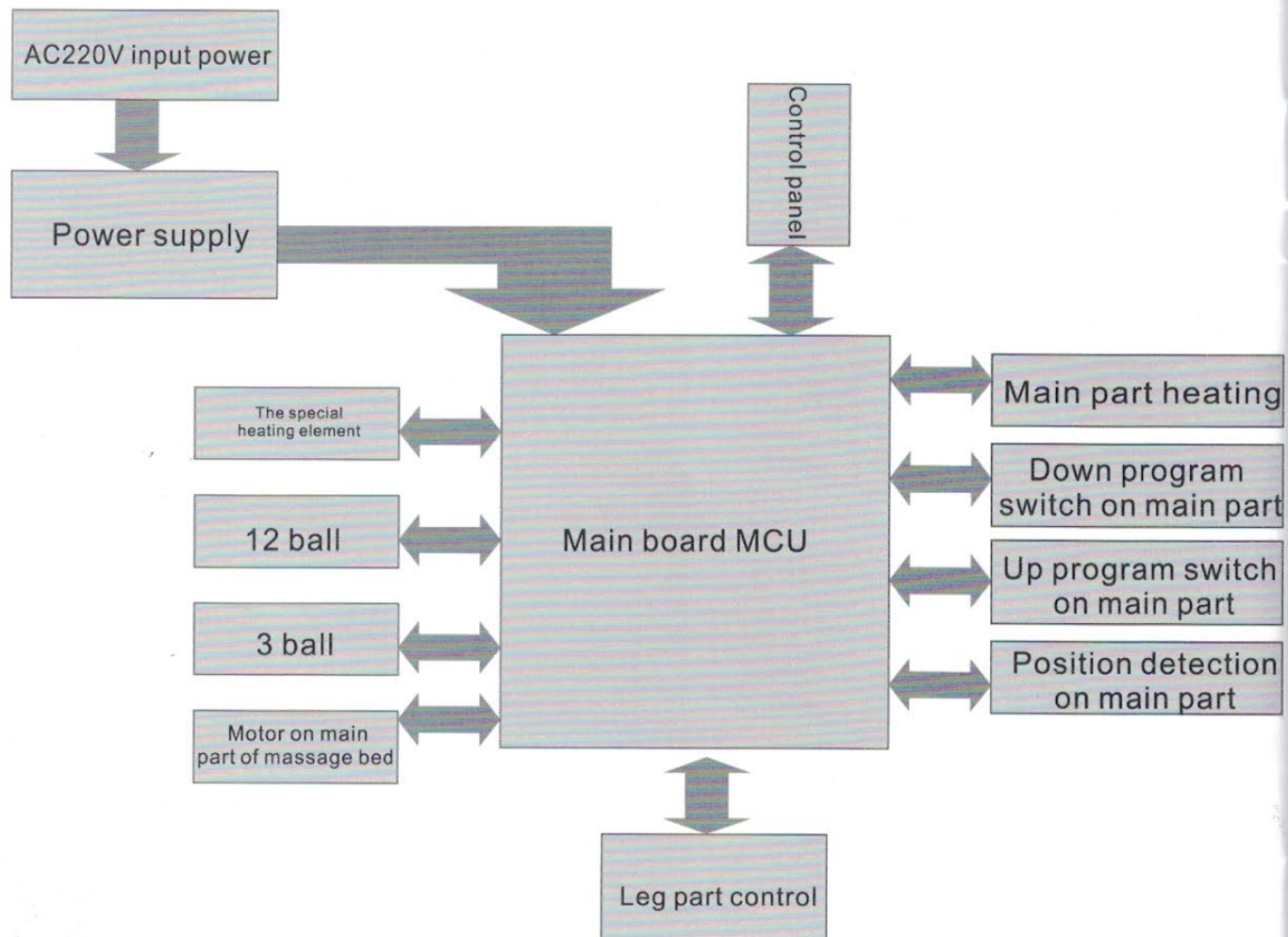

II. Phenomena of common failures and solutions

⚠ Internal machine with high voltage, the belows solutions only for professional staff's reference, non-professional staff don't disassembling machine

Phenomena of common failures	Solutions
Indication light of LCD control panel is not bright, pressing the control panel switch no reaction	<ol style="list-style-type: none"> 1.NO input AC power; please check the electric of power cord is good and connect well; whether the power switch is on; the fuse of power plug is destroyed 2.Please check the cable of LCD control panel is connected with the bed body well 3.After check with input AC power and LCD control panel light not bright, please check the internal switch output power is normal work or not, can pull out power supply and main board connect wire, then test the output power voltage is normal work or not
Temperature display error code and no heating	<ol style="list-style-type: none"> 1.Display error code E1:Means corresponding to the group of wire and accessories have a short circuit. 2.Display error code E2:Means corresponding to the group of wire and accessories have an open circuit or sockets not connect completely. 3.Display error code E3:Means communication lines between the LCD panel and main board have problem, can check connecting line of main part and the LCD panel.
Massage head of main part or leg part keep working up and down program on rail, and there is no massage head direction display on LCD	<ol style="list-style-type: none"> 1.gear wheel close to the motor is destroyed. 2.Electronic Sensitive Switch is damaged.
Massage head of main part or leg part is derailed	<ol style="list-style-type: none"> 1.The plug of up and down program switch on main part or leg part is not connect well with PCBA socket. 2.Up and down program switch is damaged.
Leg part no heating and no work	<ol style="list-style-type: none"> 1.Please check the connect line is insert to the socket well.

Note: If can't solve the failures, please don't disassembling machine by yourselves, please contact the retailer or A/S center.

